

A MODERN
MASTERPIECE

iRoll RnR 6/10
(Digital Tube Roller)

iRoll RnR 6/10 (Digital Tube Roller)

While working with all sorts of tube rollers, a few problems seem to stand out which are common among almost all devices - the use of low-quality materials and components to cut production costs results in devices that are a let down in terms of their durability and longevity. Also, due to apparent design flaws, most of these devices are large in size which is a big no-no as the space in a laboratory is already limited and these devices occupy a large space compared to their functionality.


With adequate research into user feedback and their tube rolling requirements, Neuation's R&D department helped integrate valuable insights into the design and development process of iRoll RnR 6/10, resulting in a device which comes with a reliable and durable motor for years of maintenance-free longevity. Keeping in mind the space constraints faced by most users, the exceptional design of iRoll RnR 6/10 manages to wrap incredible performance in a compact package for a wonderful user-experience.



Product Features



Brushless DC Motor
Ultra-smooth Quiet
Operation And A
Maintenance-free Long Life.

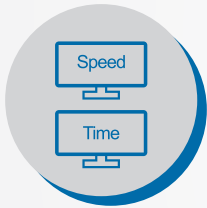

Individual Displays For
Both Time And Speed
Allows The User To Take In
All The Relevant Information
With Just One Glance.



Variable Speed Setting From
10-80 RPM With Increments
And Accuracy Of +/- 1 RPM

Variable Timer Setting
From 1 To 99 Minutes Along
With A Continuous Mode.



Allows The User To Set Upto
9 User-defined Program
For Optimum Convenience



Programmable Pulse Mode
Enables Optimum Mixing
Through Its Run-pause-run
Feature And Can Be Set
From 30 To 99 Seconds.



Product Features



Auto Power-on Feature Starts The Device Automatically After Power Failure And Runs For The Remaining Time.

Low Voltage Device Of 24v Which Makes It Safer For The User.



IP 21 Class Rating For Added Durability

Robust Construction Makes For Extremely Smooth And Stable Operation And Can Be Used In Incubators Upto 40° C And Upto 80% Humidity.



Stackable Design And Relatively Compact Footprint Saves Valuable Bench Space.

Easy Removable Roller For Convenient And Effortless Cleaning.



Circle Motion For Gentle Mixing With A Tilt Angle Of 3 Degree

Specifications

Specifications	iRoll RnR 6
Number of tubes	Upto 15 x 7 ml Blood Tubes
Speed Range	10 to 80 RPM at every 1 RPM
Motor Type	Brushless DC motor
Motion	Rocking & Rolling
Timer Range	1 to 99 mins & infinite mode
Operation Type	Continuous & Timed Operation
Removable Rollers	Yes, Tool free removal & fitment
Pulse Mode	Yes (Run-Pause-Run Mode)
Tilt Angle	+/-3°
Protection Class	IP 21
Load capacity	2 kg
Dimensions (WxDxH)	505 x 293 x 90 mm
Power Consumption	6 W
Weight	5.5 kg

Specifications	iRoll RnR 10
Number of tubes	Upto 27 x 7 ml Blood Tubes
Speed Range	10 to 80 RPM at every 1 RPM
Motor Type	Brushless DC motor
Motion	Rocking & Rolling
Timer Range	1 to 99 mins & infinite mode
Operation Type	Continuous & Timed Operation
Removable Rollers	Yes, Tool free removal & fitment
Pulse Mode	Yes (Run-Pause-Run Mode)
Tilt Angle	+/-3°
Protection Class	IP 21
Load capacity	2 kg
Dimensions (WxDxH)	502 x 448 x 93 mm
Power Consumption	6 W
Weight	7.5 kg

Ordering Information

iRoll RnR 6		
Product Name	Cat #	Details
iRoll RnR 6	TR01-1100 EU/US/UK/AU	6 Tube Roller Max RPM 80

Accessories		
Part Name	Cat #	Details
Roller for Rocking	TR01 - R1	-
Power Adaptor	TR01 - A1	-

iRoll RnR 10		
Product Name	Cat #	Details
iRoll RnR 10	TR01-2100 EU/US/UK/AU	10 Tube Roller Max RPM 80

Accessories		
Part Name	Cat #	Details
Roller for Rocking	TR01 - R1	-
Power Adaptor	TR01 - A1	-



AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

* OEM request can be entertained. Specifications can be changed without notice for quality improvement.